Sunbeam®

Renue™ Tension Relieving Heat Therapy Pad

Model 885-CN

Congratulations! You have purchased Renue™ Tension Relieving Heat Therapy Pad. Sit back, relax and start enjoying the soothing benefits of Renue™ Heat Therapy Pad.

IMPORTANT SAFETY INSTRUCTIONS - DO NOT DESTROY

WARNING: THIS PRODUCT CONTAINS MAGNETS. DO NOT USE IF YOU HAVE A PACEMAKER, ICD (implantable cardioverter defibrillator) OR AN INTERNAL INSULIN PUMP.

DANGER: TO REDUCE RISKS OF BURNS, ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

1. READ ALL INSTRUCTIONS.
2. DO NOT USE WHILE SLEEPING.
3. DO NOT USE ON AN INFANT OR ON AN ANIMAL.
4. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, A HELPLESS, SLEEPING OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, OR A PERSON WITH POOR BLOOD CIRCULATION.
5. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
6. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING. CHECK SKIN UNDER PAD FREQUENTLY.
7. DO NOT USE IN AN OXYGEN ATMOSPHERE.
8. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
9. DO NOT BUNCH, SIT ON OR AGAINST OR CRUSH PAD – AVOID SHARP FOLDS.
10. NEVER PULL THIS PRODUCT BY THE SUPPLY CORD.
11. DO NOT USE THE CORD AS A HANDLE.
12. THIS UNIT SHOULD NOT BE USED BY CHILDREN WITHOUT ADULT SUPERVISION.
13. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF THE BODY NEEDING HEAT.
14. USE THIS PAD ONLY ON A 110-120 VOLT AC CIRCUIT. UNPLUG WHEN NOT IN USE.
15. DO NOT TAMPER WITH THIS PAD IN ANY WAY.
16. THERE ARE NO USER-SERVICEABLE PARTS. IF FOR ANY REASON THIS PAD DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY FOR CONSUMER SERVICE INFORMATION.
17. DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS, ESPECIALLY ONES THAT CONTAIN HEAT-PRODUCING INGREDIENTS. SKIN BURNS COULD RESULT.
18. DO NOT BEND OR PINCH CORD AND DO NOT USE IF CORD IS DAMAGED.
19. NEVER LEAVE THE APPLIANCE UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
20. LOOP CORD LOOSELY WHEN STORING. TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS.
21. WHEN NOT IN USE, STORE THE HEATING PAD IN THE ORIGINAL PACKAGING, IN A DARK, COOL, DRY PLACE.
22. SAVE THESE INSTRUCTIONS.

This product has a polarized plug (one blade is wider than the other). This plug fits a polarized outlet only one way; if the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician.

DO NOT ATTEMPT TO DEFEAT THIS SAFETY FEATURE.

Renue™ Tension Relieving Heat Therapy is specifically designed to fit the neck and shoulder area, where users experience the most physical signs of tension and stress. Drape Renue™ Heat Therapy Pad comfortably on your shoulders making sure that the neck flap is properly in place. Weighted edges and integrated magnets help to secure Renue™ Heat Therapy Pad in place.


NOTE: This product utilizes magnets to fasten at the neck area. Remove necklaces and other metal accessories that could become attached to these magnets.

HOW TO USE YOUR RENUE™ TENSION RELIEVING HEAT THERAPY PAD

Make sure Control Connector is securely connected to the Renue™ Connector (Figure 1). Push both parts together for a tight fit. Note: The connectors are keyed to only fit one way; align arrows.

1. Plug into 120V AC electrical outlet.
2. To turn on, press Push-Button. Renue™ will NOT turn on until you press the Push-Button. **NOTE:** If the heat level indicators do not light up or stay on, the plug or Control Connector may not be firmly and completely connected. Unplug from wall outlet, check the Renue™ connection again and re-plug into outlet.
3. Select from 4 heat settings by pressing Push-Button successively to achieve desired heat level. **To turn OFF, you must scroll through the remaining heat settings by pressing the Push-Button sequentially.**
4. Allow Renue™ Heat Therapy Pad to preheat for several minutes. Renue™ has been designed to meet the maximum heat output according to UL Standards.

NOTE: A red blinking "High" light indicates an error has occurred. Unplug from outlet, verify physical condition of cord, function of AC outlet and proper connection between connector and Renue™ Pad. If blinking light persists, discontinue use and contact Consumer Services.

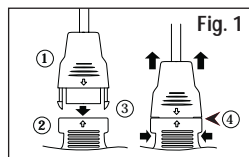1. Control Connector
2. Renue™ Connector
3. Locking Clips
4. Alignment Arrows


AUTO-OFF FEATURE

Your Renue™ controller includes an Auto-Off feature. This feature is designed to conserve energy by automatically turning off after 2 hours of continuous use. The red Auto-Off light will blink, indicating that this feature has activated. **NOTE:** After a power outage, the Auto-Off feature automatically resets the timer.

USING MOIST HEAT THERAPY

In some cases, you may find that moist heat therapy provides more comfort and tension relief. **When using moist heat therapy, you should protect furniture and clothing from becoming damp.**

1. Dampen one side of the Renue™ Heat Therapy Pad with water by patting the surface with a wet sponge or cloth or use a spray bottle to apply a fine mist. Do not thoroughly soak Renue™ Heat Therapy Pad.
2. Do not spray or dampen the connector. (Figure 2)
DANGER: Do not immerse Renue™ Heat Therapy Pad or connector in water or place under running water. Failure to comply may result in product failure and/or personal injury. **NOTE:** The electrical connector must be thoroughly dry before establishing electrical connection.

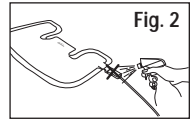WASHING INSTRUCTIONS

DO NOT DRY CLEAN! DO NOT USE BLEACH! DO NOT USE WRINGER FOR DETACHABLE PAD!

1. Disconnect pad from outlet.
2. Remove cover from pad or detach Pad Connector from Control Connector.
3. Wash by hand or as follows in automatic washer using slow agitation and normal spin speed: Fill washer with warm water. Add minimum amount of all-purpose detergent, agitate to dissolve, measure heating pad before washing. Add heating pad and soak five minutes. Agitate two minutes. Spin. Fill for cool rinse, agitate one minute then spin. Dry in automatic dryer on high heat, five minutes only, or line dry. Remove damp heating pad from dryer, block to original size, and drape to finish drying. Do not iron your heating pad.
4. Do not plug in pad and pad connector until they are completely dry.

Trouble-Shooting:

Problem	Solutions
<ul style="list-style-type: none"> • Control Blinks • No Heat • No Power 	<ul style="list-style-type: none"> • Unplug unit from wall outlet. • Verify no physical damage to the cord, controller or pad. • Verify tight connection between pad and controller. • Verify AC outlet is functional. • If blinking or failure continues, discontinue use. See Warranty Information to contact Consumer Services. • Make sure pad is not folded or bunched up in any area.
<ul style="list-style-type: none"> • Control Shuts Off 	<ul style="list-style-type: none"> • Auto Off feature is set to turn unit off after two hours of use. • Re-start unit by pressing Push-Button. • Verify tight connection between pad and controller.
<ul style="list-style-type: none"> • Pad has a Hot Spot 	<ul style="list-style-type: none"> • Make sure pad is not folded or bunched up in any area.
<ul style="list-style-type: none"> • Pad doesn't get hot enough 	<ul style="list-style-type: none"> • Re-set Heat Selection by pressing Push-Button. • Verify tight connection between pad and controller. • Heating Pad is designed to meet maximum heat output according to UL Standards.

Note: this heating pad has a built in safety feature that prevents the heating pad from being turned on immediately after the auto off feature has turned it off. Please wait at least 1 minute before turning your heating pad back on.

5-YEAR LIMITED WARRANTY

Sunbeam Products, Inc. doing business as Jarden Consumer Solutions or if in Canada, Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions (collectively "JCS") warrants that for a period of five years from the date of purchase, this product will be free from defects in material and workmanship. JCS, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty. This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. JCS dealers, service centers, or retail stores selling JCS products do not have the right to alter, modify or in any way change the terms and conditions of this warranty.

This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use of improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than JCS or an authorized JCS service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

What are the limits on JCS's Liability?

JCS shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty. JCS disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. JCS shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party. Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How to Obtain Warranty Service

In the U.S.A.

If you have any question regarding this warranty or would like to obtain warranty service, please call 1 800 435-1250 and a convenient service center address will be provided to you.

In Canada

If you have any question regarding this warranty or would like to obtain warranty service, please call 1 800 667-8623 and a convenient service center address will be provided to you. In the U.S.A., this warranty is offered by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions located in Boca Raton, Florida 33431. In Canada, this warranty is offered by Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions, located at 20B Hereford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE
©2010 Sunbeam Products, Inc. doing business as Jarden Consumer Solutions. All rights reserved. Distributed by Sunbeam Corporation (Canada) Limited, doing business as Jarden Consumer Solutions, 20B Hereford Street, Brampton, Ontario L6Y 0M1.

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