Trouble-Shooting:

Problem	Solutions
<ul style="list-style-type: none">No Heat No Power	<ul style="list-style-type: none">Unplug unit from wall outlet. Verify no physical damage to the cord, controler or pad. Verify tight connection. Verify AC outlet is functional. Restart unit by sliding control to off position then back to desired heat level. Make sure pad is not folded or bunched in any area
<ul style="list-style-type: none">Control Shuts Off	<ul style="list-style-type: none">Auto Off feature is set to turn unit off after two hours of use. Re-start unit by sliding control to off positio then back to desired heat level Verify tight connection
<ul style="list-style-type: none">Pad has a Hot Spot Pad doesn't get hot enough	<ul style="list-style-type: none">Make sure pad is not folded or bunched up in any area Re-set Heat Selection by pressing Push-Button Verify light connection. Heating Pad is designed to meet maximum heat output according to UL Standards.

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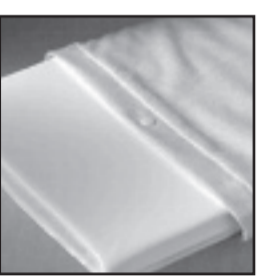
Figure 1

Figure 2

FABRIC COVER WASHING INSTRUCTIONS

- Disconnect pad from outlet.
- Remove cover from pad.
- Machine or hand wash; use mild detergent. Wash with cool water only.
- Hang to dry or machine dry on gentle, cool cycle. (Figure 2)

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7-Year Limited Warranty

Sunbeam Products, Inc. or if in Canada, Sunbeam Corporation (Canada) Limited (collectively "Sunbeam") warrants that for a period of seven years from the date of purchase, this product will be free from defects in material and workmanship. Sunbeam, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty.

This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. Sunbeam dealers, service centers, or retail stores selling Sunbeam products do not have the right to alter, modify or any way change the terms and conditions of this warranty.

This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than Sunbeam or an authorized Sunbeam service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

What are the limits on Sunbeam's Liability?

Sunbeam shall not be liable for any incidental or consequential damages caused by the breach of any expressed, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty.

Sunbeam disclaims all other warranties, conditions or representations, expressed, implied, statutory or otherwise.

Sunbeam shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamen-tal or otherwise, or for any claim brought against purchaser by any other party.

Some provinces, states or jurisdictions do not allow the exclusion or limitation of inci-dental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How to Obtain Warranty Service

If you have any questions regarding this warranty or would like to obtain warranty service, please call and a convenient service center address will be provided to you, or visit our web site at www.sunbeam.com

In the U.S.A. - This warranty is offered by Sunbeam Products, Inc., located in Boca Raton, Florida 33431. 1-800-435-1250

In Canada, this warranty is offered by Sunbeam Corporation (Canada) Limited, located at 5975 Falbourn Street, Mississauga, Ontario L5R 3V8 1-800-667-8623

©2007 Sunbeam Products, Inc. doing business as Jarden Consumer Solutions. All rights reserved.
Printed in China

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IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

DANGER: TO REDUCE RISKS OF BURNS, ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- READ ALL INSTRUCTIONS
- DO NOT USE WHILE SLEEPING.
- DO NOT USE ON AN INFANT OR ON AN ANIMAL.
- THIS PAD IS NOT TO BE USED ON OR BY AN INVALID,A SLEEPING OR UNCONSCIOUS PERSON,A PERSON WITH DIABETES,OR A PERSON WITH POOR BLOOD CIRCULATION.
- DO NOT USE ON AREAS OF INSENSITIVE SKIN.
- BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING.CHECK SKIN UNDER PAD FREQUENTLY.
- DO NOT USE IN OXYGEN ENRICHED ENVIRONMENT.
- DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
- DO NOT SIT ON OR AGAINST OR CRUSH PAD-AVOID SHARP FOLDS.
- NEVER PULL THIS PAD BY THE SUPPLY CORD
- DO NOT USE TEH CORD AS A HANDLE.
- CAREFULLY EXAMINE INNER COVER BEFORE EACH USE.DISCARD THE PAD IF INNER COVERING SHOWS ANY SIGNS OF DETERIORATION.
- THIS UNIT SHOULD NOT BE USED BY CHILDREN WITHOUT ADULT SUPERVISION.
- PLACE PAD ON TOP OF AND NOT UNDER THE PART OF THE BODY NEEDING HEAT.
- USE THIS PAD ONLY ON A 110-120 VOLT AC CIRCUIT. UNPLUG WHEN NOT IN USE.
- DO NOT TAMPER WITH THIS PAD IN ANY WAY.
- THERE ARE NO USER-SERVICEABLE PARTS.IF FOR ANY REASON THIS PAD DOES NOT FUNCTION SATISFACTORILY,SEE WARRANTY FOR CONSUMER SERVICE INFORMATION.
- DO NOT USE THIS PAD WITH LINIMENT,SALVE OR OINTMENT PREPARATIONS,ESPECIALLY ONES THAT CONTAIN HEAT-PRODUCING INGREDIENTS.SKIN BURNS COULD RESULT.
- DO NOT BEND OR PINCH CORD AND DO NOT USE IF CORD IS DAMAGED.
- NEVER LEAVE THE APPLIANCE UNATTENDED,ESPECIALLY IF CHILDREN ARE PRESENT.
- SAVE THESE INSTRUCTIONS.

FIGURE 1	FIGURE 2
-----------------	-----------------

This appliance has a polarized plug (one blade is wider than the other).This plug fits a polarized outlet only one way; if the plug does not fit fully into the outlet, reverse the plug.If it still does not fit,contact a qualified electrician. DO NOT ATTEMPT TO DEFEAT THIS SAFETY FEATURE.

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
Model / Modèle 756-505

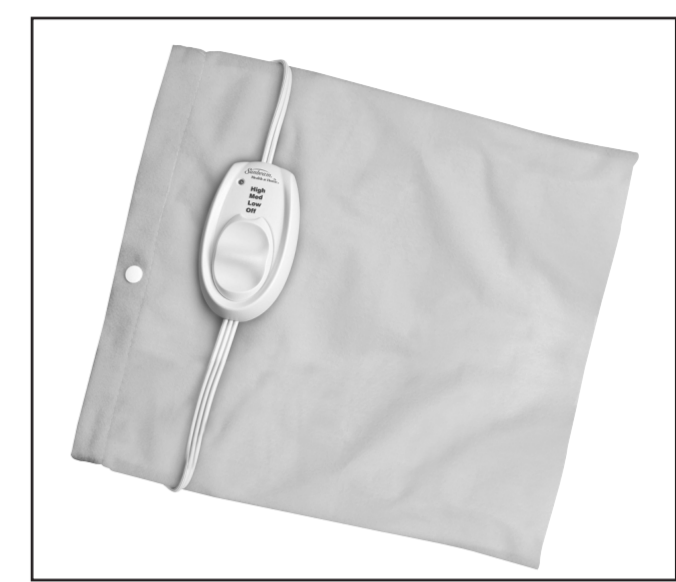IMPORTANT NOTES

Please read this entire manual carefully before using this device. Product may vary from picture.

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FOR MODELS WITH GEL PACK:

- Heat gel in microwave.
 - Lay flat in microwave
 - Make sure gel is spread evenly across pack for proper heating.
- CAUTION: UNEVEN GEL COULD CAUSE HOT SPOTS**
 - Heat at full power (700 Watts) for one minute.
 - Remove and check heat guide.
- If not hot enough, return to microwave and heat in ten second intervals, checking Heat Guide until ready.
- If Heat Guide reads TOO HOT, let cool before use.

Gel pack safety reminders

- Gel Pack can be used without heating pad, but must be enclosed in cover.
- NEVER HEAT PAD FOR MORE THAN 90 SECONDS AT ONE TIME OR GEL PACK MAY RUPTURE AND LEAK**
- NEVER BOIL GEL PACK**
- insert Heating Pack and Gel Pack in cover
 - Lay heating pad flat.
 - Attach Gel Pack to heating pad. NOTE:You can use heating pad WITHOUT Gel pack.
 - Slip inside cover, making sure cord and controls remain on outside.
 - Close snap
- Plug cord in wall.

NOTE: Although your heating pad is meant to get hot, it should NEVER hurt or feel like it is burning your skin. If it does, remove the heating pad immediately.

For models with auto off feature:

Heating pad will automatically shut itself off after 2 hours of use. Restar unit by sliding to off position then back to desired heat level.

HOW TO USE YOUR HEATING PAD:

- Lay the cloth cover on a flat surface with the printing face up. Open the snap at the bottom.
- Insert the pad into the cover completely, fitting one side at a time. Make sure the cord and control remain on the outside of the cover.
- The pad will lie flat, filling the entire cover. Close snap.
- Plug into any 120 volt AC electrical outlet and set control to the desired heat.

WARNING: Never use the heating pad without the cloth cover.

USING MOIST HEAT OPTION ON SELECTED MODELS
Your heating pad may include a separate sponge for moist heat.

- Run warm water over sponge until it is thoroughly saturated.
- Gently wring out excess water.
- Place the sponge on top of the heating pad and slip the cloth cover over the heating pad. Snap cover closed.

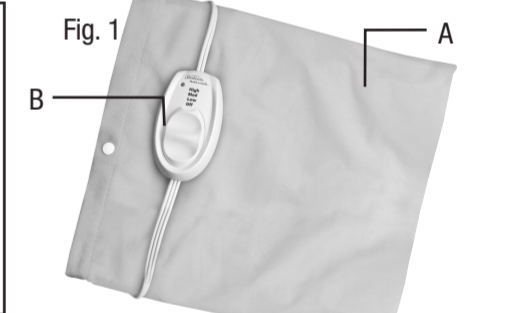
FIGURE 1	FIGURE 2
-----------------	-----------------

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