Sunbeam

Renue™ Massaging
Heating Pad

Models: 730-811-CN
730-805-CN
730-811A-CN

IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

DANGER: TO REDUCE RISKS OF BURNS, ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:


1. READ ALL INSTRUCTIONS.
2. DO NOT USE WHILE SLEEPING.
3. DO NOT USE ON AN INFANT OR ON AN ANIMAL.
4. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, A SLEEPING OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, OR A PERSON WITH POOR BLOOD CIRCULATION.
5. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
6. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING. CHECK SKIN UNDER PAD FREQUENTLY.
7. NEVER USE PAD WITHOUT COVER IN PLACE.
8. DO NOT USE IN OXYGEN ENRICHED ENVIRONMENT.
9. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
10. DO NOT SIT ON OR AGAINST OR CRUSH PAD; AVOID SHARP FOLDS.
11. NEVER PULL THIS PAD BY THE SUPPLY CORD.
12. DO NOT USE THE CORD AS A HANDLE.
13. CAREFULLY EXAMINE INNER COVER BEFORE EACH USE. DISCARD THE PAD IF INNER COVERING SHOWS ANY SIGNS OF DETERIORATIONS.
14. THIS UNIT SHOULD NOT BE USED BY CHILDREN WITHOUT ADULT SUPERVISIONS.
15. USE THIS PAD ONLY ON A 110-120 VOLT AC CIRCUIT. UNPLUG WHEN NOT IN USE. DO NOT USE WITH GENERATORS, POWER CONVERTERS, OR INVERTERS.
16. DO NOT TAMPER WITH THIS PAD IN ANY WAY.
17. THERE ARE NO USER-SERVICEABLE PARTS. IF FOR ANY REASON THIS PAD DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY FOR CONSUMER SERVICE INFORMATION.
18. DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS, ESPECIALLY ONES THAT CONTAIN HEAT PRODUCING INGREDIENTS, SKIN BURNS COULD RESULT.
19. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF BODY NEEDING HEAT.
20. LOOP CORD LOOSELY WHEN STORING, TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS.
21. SAVE THESE INSTRUCTIONS.

WARNING: Do not massage abdomen when pain is present. Do not massage areas that are swollen, inflamed, or covered with skin eruptions. Do not use on unexplained calf pain. Consult a physician before using any equipment of this type. INCLUDES A UL/CUL LISTED TRANSFORMER ELECTRIC HEATING PAD – FOR HOUSEHOLD USE ONLY

THIS APPLIANCE HAS A POLARIZED PLUG (ONE BLADE IS WIDER THAN THE OTHER). TO REDUCE THE RISK OF ELECTRIC SHOCK, THIS PLUG IS INTENDED TO FIT IN A POLARIZED OUTLET ONLY ONE WAY. IF THE PLUG DOES NOT FIT FULLY IN THE OUTLET, REVERSE THE PLUG. IF IT STILL DOES NOT FIT, CONTACT A QUALIFIED ELECTRICIAN TO INSTALL THE PROPER OUTLET.

DO NOT ATTEMPT TO MODIFY THIS PLUG OR DEFEAT THIS SAFETY FEATURE IN ANY WAY.

- A. Soft Cover
- B. Heat Control Lo/Hi
- C. Massage/Vibration Intensity Lo/Hi
- D. Power Adapter Plug Port
- E. Power Adapter

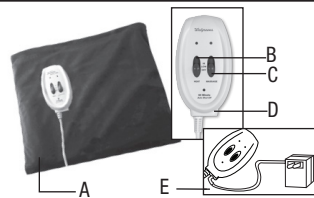FEATURES

Your Sunbeam® Renue™ Massaging Heating Pad comes with the following:

- Soft cover that is removable and machine washable.
- Slide control with 2 heat and 2 massage settings.

Features Of Your Renue™ Massaging Heating Pad

- Soothing warmth and vibration for sore muscles and temporary arthritic pain
- Use heat and massage therapy independently or together
- 2 heat and massage levels
- Moist heat option
- Slide control with ON indicator light
- Auto shut-off after one hour

How to Use the Renue™ Massaging Heating Pad

For Dry Heat Application

1. Lay the cloth cover on a flat surface with the embroidered side up.
2. Insert the pad, with the heated side up, into the cover completely, fitting one side at a time.
3. Make sure the cord and controller are on left side, and both remain on the outside of the cover.
4. The pad will lie flat, filling the entire cover.
5. Plug into nearest power supply.
6. Plug power adapter into controller at plug port (D).
7. Turn the controller on to desired heat and/or massage setting. Allow 15 minutes for the pad to reach the high heat setting.
8. Relax and enjoy. Do not use for more than 60 minutes at a time (not counting warm-up time).
9. When finished, turn controller off and unplug the heating pad.

CAUTION: NEVER USE PAD WITHOUT COVER IN PLACE.

Fabric Cover Washing Instructions

CAUTION: DO NOT DRY CLEAN! DO NOT BLEACH! DO NOT USE WRINGER!

1. Disconnect pad from outlet.
2. Remove cover from pad.
3. Machine wash cover in cool water on gentle cycle, or hand wash. Use mild detergent.
4. Hang to dry or machine dry on gentle, cool cycle.

2-YEAR LIMITED WARRANTY

Sunbeam Products, Inc. doing business as Jarden Consumer Solutions or if in Canada, Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions (collectively "JCS") warrants that for a period of 2 years from the date of purchase, this product will be free from defects in material and workmanship. JCS, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty. This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. JCS dealers, service centers, or retail stores selling JCS products do not have the right to alter, modify or in any way change the terms and conditions of this warranty. This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use of improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than JCS or an authorized JCS service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

What are the limits on JCS's Liability?

JCS shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty.

JCS disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. JCS shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party. Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How to Obtain Warranty Service

In the U.S.A. - If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-435-1250 and a convenient service center address will be provided to you.

In Canada, - If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-667-8623 and a convenient service center address will be provided to you.

In the U.S.A., this warranty is offered by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions located in Boca Raton, Florida 33431. In Canada, this warranty is offered by Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions, located at 20 B Herford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.

© 2016 Sunbeam Products, Inc. doing business as Jarden Consumer Solutions. All rights reserved. In the U.S., distributed by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions, Boca Raton, FL 33431. In Canada, imported and distributed by Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions, 20 B Herford Street, Brampton, Ontario L6Y 0M1. Sunbeam Products, Inc. is a subsidiary of Newell Brands Inc. (NYSE: NWL).

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